


Russell PneumoFix Thoracic Decompression Needle (12 ga x 4.33")



The Russell PneumoFix® is a sterile chest decompression device designed for the management of tension pneumothorax, simple pneumothorax and pleural effusion. Designed by practicing clinicians who understand the limitations of conventional equipment, the Russell PneumoFix® has many features that make it a truly useful tool for the emergency care provider.

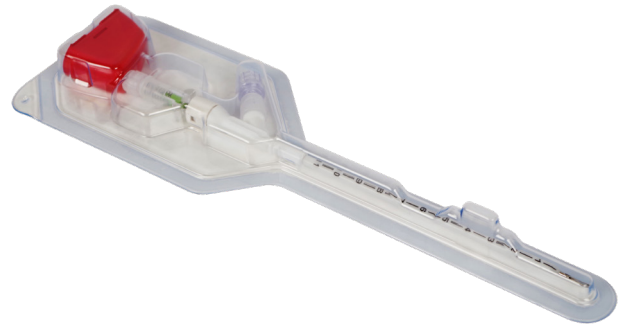

Key Benefits:

- Simple to use, minimal risk involved.
- Designed for efficient insertion to be inserted without the use of a scalpel or skin incision.
- Veress-tipped needle to reduce the risk of injury to lung tissue.
- 11cm long catheter long enough to reach the pleural cavity of the vast majority of patients.
- Low pressure one-way valve to allow air and fluid out but nothing in.
- Graduated markings and X-ray detectable catheter for depth recording and accurate localization at hospital.
- 12-Gauge catheter allows for rapid venting of air in tension pneumothorax.

Specifications

- PN/SKU: PDF112
- Product code: PDF102
- Pack Dimensions: 28.00cm x 8.50cm
- Product shelf life: 5 years

For further information please contact us at:

Tel: 855-428-6074 Email: info@safeguardmedical.com

Visit: safeguardmedical.com